

ONE SOURCE. INFINITE SOLUTIONS.®

Vacuum Pumps



Part of Thermo Fisher Scientific

For filtration, degassing, desiccation, distillation, concentration, drying...
Precisely tuned to your laboratory applications.



Fisher Scientific presents:

Dr Blue



Fisher Scientific Vacuum Pumps

Pump Selection and Application Guide

Perfectly tuned for

- Filtration
- Degassing
- Desiccation
- Distillation
- Concentration
- Drying ...

Innovative Technology: Fisher Scientific vacuum pumps combine high tech materials with outstanding engineering to provide:

- Long lifetime
- Low consumption
- Small footprint
- Reduced noise level

For Multiple Apps: Fisher Scientific vacuum pumps are designed to fit most laboratory applications, both routine and research.

- High pumping speed
- Wide ultimate vacuum range
- Excellent chemical resistance
- Minimum maintenance

Quick Vacuum Pump Selection and Application Guide

General application		FB65540	FB65451	FB65452	FB65453	FB65454	FB65458	FB65459	FB65460	FB65461	FB65464	see page
Selecting factor	Ultimate vacuum											
Operating vacuum range	Pressure above atmospheric	X										3
	Vacuum to 290 mbar	X										3
	Vacuum to 60 mbar		X									4
	Vacuum to 8 mbar			X	X							4
	Vacuum to 2 mbar					X						4
	Vacuum to 10 ⁻² mbar						X					5
	Vacuum to 10 ⁻³ mbar							X	X	X	X	5



Dr Blue says...

This guide covers most standard applications for vacuum pumps. However, if your application is not mentioned, we may still have a pump to meet your requirement. Please contact your local Fisher Scientific representative or Technical Support.

Pump Selection / Application Chart

Application		FB65540	FB65451	FB65452	FB65453	FB65454	FB65458	FB65459	FB65460	FB65461	FB65464	see page
Filtration												
Vacuum or pressure method	Pressure method	X										3
	Vacuum method	X	X									3
Gas transfer												
Low volume		X	X									3
Desiccators and pistol drying												
Residual moisture level	Partial drying	X										3
	Full drying		X	X								3
Rotary evaporator												
Boiling point of solvent used	Low boiling point			X	X							4
	High boiling point (DMF)					X						4
Centrifugal concentrator												
Boiling point of solvent used	Low boiling point				X							4
	High boiling point (DMF)					X						4
Vacuum oven												
Oven chamber volume	Around 25 liters				X							4
	Around 55 liters				X	X						4
Freeze dryer												
	up to 2 kg							X				5
	up to 4 to 6 kg								X			5
Ice holding capacity	up to 8 kg									X		5
	up to 16 kg										X	5
Glove box												
Standard for all sizes										X		5
High vacuum distillation												
Standard for all sizes							X	X				5

Fisher Scientific Vacuum Pumps

General Duty Fisher Scientific Vacuum Pumps

Fisher Scientific oil-free piston pumps and diaphragm pumps for a variety of applications such as filtration of low volumes, degassing, gas transfer and sampling.

Piston Pump - FB65540

The FB65540 piston pump is ideal for the basic filtration of non-aggressive solvents.

Features:

- Lightweight and compact
- Covers most standard applications
- Economically priced but with high performance
- Quick and simple to use due to stylish design
- Portable
- Certifications: CE marked, UL/CSA approved



Diaphragm Pump - FB65451

The FB65451 double head diaphragm pump offers a choice of two speeds for softer filtration and is also ideal for drying and gas transfer.

Features:

- Two speed performance for a gentle vacuum
- Compact, lightweight and portable
- Quiet operation and very low vibration
- Use in lab or field
- Certifications: CE marked, UL/CSA approved



Technical Features:

Type	FB65540	FB65451
Ultimate vacuum	< 290 mbar	< 60 mbar
Pumping speed at 50/60Hz	9.2 l/min	10 l/min
max. overpressure	2.3 bar	0.3 bar
Inlet and outlet connection	hose nozzle DN 6	hose nozzle DN 6/8, silencer
Sound level	< 56 db (A)	< 42 db (A)
Motor power	max. 160 W	20 W
Type of protection (motor)	IP 48	IP 44
Weight	2.3 kg	2.3 kg
Overall dimensions	191/102/183 mm	115/165/145 mm

Ordering Information:

Type	Cat. No.
Piston pump 230 V 50 Hz	FB65540
Piston pump 115 V 50/60 Hz	FB65540A
Diaphragm pump 90-240 V 50/60 Hz, incl. AC/DC adapter	FB65451
Maintenance kit	402031

Dr Blue says...

For a convenient filtration process where you always have the right vacuum, the FB65451 is ideal. FB65540 and FB65451 are not recommended for acid and strong chemical vapors. For a chemical resistant model please see page 4.



Fisher Scientific Vacuum Pumps


Chemical Resistant Diaphragm Vacuum Pumps

This range of Fisher Scientific diaphragm pumps is chemically resistant, making it ideal for applications involving aggressive solvents and acidic vapors. They are an excellent replacement for rotary vane pumps and water aspirators, especially in rough vacuum applications where corrosive vapors are encountered. The diaphragms and all gas-contaminated parts are made from PTFE and PTFE compound.

The pumps and connection heads are carbon fiber reinforced to provide electrical conductivity and prevent electrostatic charging. A long lifetime and minimal maintenance make these pumps a great choice. Other features include an integrated protective switchgear for the motor, hose nozzles for 8 mm ID tubing and vibration damping feet.

Features:

- Handle corrosive and aggressive vapour without damage
- Clean, reliable vacuum pumps to suit your needs
- Innovative design for excellent running
- Compact and easy to use
- Certifications: CE marked, UL/CSA approved

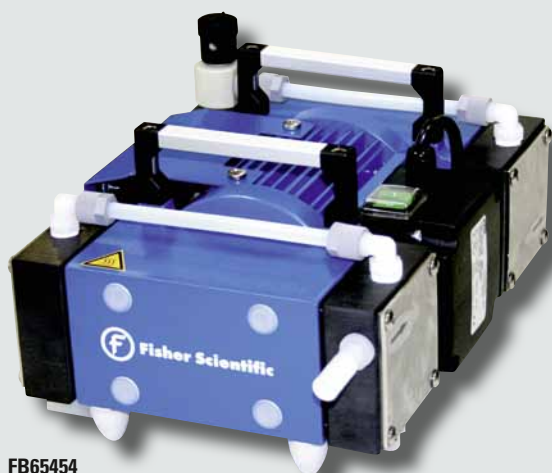

FB65454



FB65453



FB65452

Technical Features:

Type	FB65452	FB65453	FB65454
Number of heads / stages	2/2	2/2	4/3
Ultimate vacuum	< 8 mbar	< 8 mbar	< 2 mbar
Pumping speed at 50/60Hz	16/18 l/min	38/41 l/min	33/36 l/min
Max. overpressure	1 bar	1 bar	1 bar
Inlet and outlet connection	Hose nozzle DN 8	Hose nozzle DN 8	Hose nozzle DN 8
Sound level	< 44 db (A)	< 44 db (A)	< 44 db (A)
Motor power	60W	180 W	90 W
Type of protection (motor)	IP 54	IP 54	IP 54
Weight	6.5 kg	11.2 kg	10.3 kg
Overall dimensions (W/D/H)	195/235/145 mm	230/265/169 mm	200/260/150 mm



Dr Blue says...

Check the pump on a regular basis for unusual running noises and heat building up on the surface of the pump. We recommend replacing parts using the maintenance kit after 10,000 operating hours or less, depending upon the application.

Ordering Information:

For accessories see page 6.

Type	Cat. No.	Cat. No.	Cat. No.
Diaphragm Vacuum Pump 230 V 50/60 Hz	FB65452	FB65453	FB65454
Diaphragm Vacuum Pump 115 V 50/60 Hz	FB65452A	FB65453A	FB65454A
Maintenance kit	402008	402041	402015

Fisher Scientific Vacuum Pumps

Two Stage Rotary Vane Vacuum Pumps

These two-stage rotary vane pumps provide a high water vapor tolerance. They are designed for a wide range of laboratory applications and continuous running. The low weight, the small footprint and the high pumping speed, as well as the increased chemical resistance are advantageous features in research and routine processes.

Features:

- New compact design, low weight and easy to use
- Quiet running and excellent ultimate vacuum
- Active corrosion protection
- Easy maintenance
- Certifications: CE marked, UL/CSA approved



FB65459 / FB65460



FB65458



FB65461 / FB65464



Dr Blue says...

All Fisher Scientific rotary vane pumps are supplied with oilfilling of Labovac 10 oil. This enables best performance from your pump. Labovac 10 oil has a low vapor pressure and is perfect for general use.

Other oils are available for more special applications. Please contact your local Fisher Scientific representative or Technical Support.

Technical Features:

Type	FB65458	FB65459	FB65460	FB65461	FB65464
Number of stages	2	2	2	2	2
Ultimate vacuum					
- without gas ballast total	1 x 10 ⁻² mbar	2 x 10 ⁻³ mbar	2 x 10 ⁻³ mbar	2 x 10 ⁻³ mbar	2 x 10 ⁻³ mbar
- with gas ballast total	0.5 mbar	1 x 10 ⁻² mbar	1 x 10 ⁻² mbar	1 x 10 ⁻² mbar	2 x 10 ⁻² mbar
Pumping speed at 50/60Hz	30/36 l/min	77/92 l/min	97/110 l/min	120/143 l/min	183/220 l/min
Max. allowed inlet / outlet pressure	1 bar	1 bar	1 bar	1 bar	1 bar
Inlet and outlet connection	DN 16 KF	DN 16 KF	DN 16 KF	DN 25 KF	DN 25 KF
Sound level	< 48 db (A)	< 49 db (A)	< 49 db (A)	< 49 db (A)	< 49 db (A)
Motor power	120 W		200 W		550 W
Type of protection (motor)	IP 54		IP 54		IP 54
Water vapor tolerance	33 mbar		60 mbar		35 mbar
Oil capacity	250 ml	530 ml	550 ml	600 ml	1000 ml
Weight	8.0 kg	17.5 kg	19.5 kg	21.5 kg	35.0 kg
Overall dimension (W/D/H)	300/120/150 mm	415/150/235 mm	430/150/235 mm	445/150/235 mm	540/170/260 mm

Ordering Information:

For accessories see page 6.

Type	Cat. No.	Cat. No.	Cat. No.	Cat. No.	Cat. No.
Rotary Vane Vacuum Pump 100-120 V/200-240 V 50/60 Hz		FB65459 ¹	FB65460 ¹	FB65461 ¹	FB65464 ¹
Rotary Vane Vacuum Pump 230 V 50/60 Hz	FB65458 ²				
Rotary Vane Vacuum Pump 115 V 50/60 Hz	FB65458A ²				
Maintenance kit	302012-1	302075	302076	302077	302079

Included in the delivery range:

¹ Two-stage rotary vane vacuum pump complete with anti-suction back valve, gas ballast valve, oil charging, centering ring, clamping ring and motor overload protection

² Two-stage rotary vane vacuum pump with oil charging, clamping ring, centering ring, motor protection

Fisher Scientific Vacuum Pumps

Accessories

Accessories for Chemical Resistant Diaphragm Vacuum Pumps

Fitting sets for Chemically Resistant Diaphragm Pumps FB65452, FB65453, FB65454:

- All in one kits, e.g. for connection to a vacuum oven with DN 16 KF or DN 25 KF flange. Each kit comes ready to use.

Type	Connection flange	Cat. No.
Fitting Set 1 consisting of - hose nozzle DN 16-8 - centering ring, clamping ring DN 16 KF - vacuum tube 18/8x5 (2.5 m)	DN 16 KF	404002
Fitting Set 2 consisting of - hose nozzle DN 16-8 - centering ring, clamping ring DN 25 KF - vacuum tube 18/8x5 (2.5 m)	DN 25 KF	404004



Fitting set

Accessories for Two Stage Rotary Vane Vacuum Pumps

Fitting sets for Two Stage Rotary Vane Vacuum Pumps FB65459, FB65460, FB65461, FB65464:

- All in one kits, e.g. for connection to freeze dryers with DN 16 KF or DN 25 KF flange. Each kit comes ready to use.

Type	Connection flange	Cat. No.
Fitting Set 1 consisting of - 2 way ball valve - clamping ring, centering ring DN 16 KF - PVC tubing (1.5 m) - condenser AKD 16	DN 16 KF	330001
Fitting Set 2 consisting of - 2 way ball valve - clamping ring, centering ring DN 25 KF - PVC tubing (1.5 m) - condenser AKD 25	DN 25 KF	330002



Fitting set

Exhaust filters for Two Stage Rotary Vane Vacuum Pumps FB65458, FB65459, FB65460, FB65461, FB65464:

These compact filters can be simply mount on the exhaust port to collect oil expelled during operation and to reduce pump noise.

A replaceable filter element within the oil filter separates the oil from particles to collect them in the filter housing. When the filter element becomes saturated simply unscrew the filter for easy replacement.

- Condenser AKD for effective oil mist elimination (not for FB65458)
- Oil mist filter OME for economical oil mist elimination

Type	Connection flange	Cat. No.
Condenser AKD 16	DN 16 KF	300883
Condenser AKD 25	DN 25 KF	300885
Oil mist filter OME 10/16	DN 16 KF	700010
Oil mist filter OME 10/25	DN 25 KF	700011



Condenser AKD



Oil mist filter OME



Dr Blue says...

Looking for further accessories such as inlet separators, cold traps, vacuum gauges and sensors etc?

If so, please contact your local Fisher Scientific representative or Technical Support and they will be happy to help.

1. Characteristic of the application

Specify type and/or size

Rotary evaporator

Vacuum oven

Vacuum concentrator

Freeze dryer

High Vacuum distillation

Vacuum filtration

Pressure filtration

Desiccators and pistol drying

Glove boxes

Others

Company

Name

Address

ZIP / City

Country

Phone / Fax

E-mail



Dr Blue says...

If the above pumps still do not meet your requirements, then please send us full details of your application using the Pump Questionnaire. If you have a drawing or other details, feel free to send them with the questionnaire. Alternatively, please contact your local Fisher Scientific representative or Technical Support. We will then endeavor to find a pump to meet your requirement.

2. Characteristic of the gas

Description

Description of the gas

Temperature of the gas

Volume of the gas

3. Characteristic of the pump

Description

Pumping speed

Ultimate vacuum to

Connection of the vacuum – inlet and outlet

Voltage and frequency of the motor

Protection class of the motor

4. Remarks

Australia

Thermo Fisher Scientific
5 Caribbean Drive
Scoresby
Victoria 3179 Australia
Tel: +61 3 9757 4300
Fax: +61 3 9763 2791
www.thermofisher.com.au

China

Fisher Scientific Worldwide Shanghai Co., Ltd.
20/F Longlife Mansion
No. 1566 YanAn West Road
200021 Shanghai
China
Tel: +86 21 5258 1100
Fax: +86 21 5258 0114
www.fishersci.com.cn

France

Fisher Scientific
Parc d'Innovation, BP 50111
Boulevard Sebastien Brant
67400 Illkirch-Cedex
France
Tel: +33 (0)3 88 67 14 14
Fax: +33 (0)3 88 67 11 68
www.fr.fishersci.com

Japan

Fisher Scientific Japan Ltd.
Tokyo Nissan Bldg. 3F 3-25-9
Taishi-do Setagaya-ku
154-0004 Tokyo
Japan
Tel: +81 334-196778
Fax: +81-334-196781
www.fishersci.co.jp

New Zealand

Biolab Limited
Thermo Fisher Scientific
244 Bush Road
Albany
Auckland New Zealand
Tel: +64 09 980 6700
Fax: +64 09 980 6788
www.thermofisher.co.nz

South Korea

Fisher Scientific Korea Limited
Sambu Bldg. 13F, 676
Yoksam-Dong Kangnam-Ku
140-210 Seoul, South Korea
Tel: +82 2 527 0300
Fax: +82 2 527 5662
www.fishersci.co.kr

UK

Fisher Scientific UK Limited
Bishop Meadow Road
Loughborough
Leicestershire LE11 5RG UK
Tel: +44 150 923 1166
Fax: +44 150 923 1893
www.fisher.co.uk

Austria

Fisher Scientific Austria GmbH
Rudolf von Alt-Platz 1
1030 Wien
Austria
Tel: +43 17685611
Fax: +43 17685612
www.at.fishersci.com

Czech Republic

Fisher Scientific spol. s r.o.
Kosmonautu 324
53009 Pardubice
Czech Republic
Tel: +42 (0) 46 6 43 50 01
Fax: +42 (0) 46 6 435 008
www.fishersci.cz

Germany

Fisher Scientific GmbH
Im Heiligen Feld 17
58239 Schwerte
Germany
Tel: +49 (0) 2304 932 5
Fax: +49 (0) 2304 932 950
www.de.fishersci.com

Latin America

Fisher Scientific Latin America
Headquarters
3970 John's Creek Court
Suite 500
Suwanee
GA 30024 USA
Tel: +1 770-871-4725
Fax: +1 770-871-4726
www.fishersci.com

Norway

Fisher Scientific
Frysjaeveien 33E
0884 Oslo
Norway
Tel: +47 22 95 59 43
Fax: +47 22 95 59 40
www.fishersci.no

Spain

Fisher Scientific SL
C/ Luis I 9
28031 Madrid
Spain
Tel: 91 380 67 10
Fax: 91 380 85 02
www.es.fishersci.com

Belgium

Fisher Scientific SPRL
BP 567
7500 Tournai 1
Belgium
Tel: +32 056 260 260
Fax: +32 056 260 270
www.be.fishersci.com

Denmark

Fisher Scientific
Postboks 60
Industrivej 3
3550 Slangerup
Denmark
Tel: +45 70 27 99 20
Fax: +45 70 27 99 29
www.fishersci.de

India

Fisher Scientific India
Godrej Coliseum 101A-101B
Somaiya Hospital Road
Off Eastern Express Highway, Sion (East)
IN-400 022 Mumbai
Tel: 022-6680 3000
Fax: 022-6680 3001/02
www.fishersci.in

Malaysia

Fisher Scientific (M) Sdn. Bhd.
No. 3 Jalan Sepadu 25/123
Taman Perindustrian Axis Seksyen 25
40400 Shah Alam Selangor Darul Ehsan
Malaysia
Tel: +60 3-5122-8888
Fax: +60 3-5121-8899
www.fishersci.com.my

Portugal

Fisher Scientific
Rua Pedro Álvares Cabral nº24
3ºD Edifício Euro
Infantado 2670-391 Loures
Portugal
Tel: +351 21 425 33 50/4
Fax: +351 21 425 33 51
www.pt.fishersci.com

Sweden

Fisher Scientific
Box 9193
400 94 Göteborg
Sweden
Tel: +46 31 352 32 00
Fax: +46 31 352 32 50
www.fishersci.se

Canada

Fisher Scientific Limited
Ottawa Admin. Center
112 Colonnade Road
Nepean
Ontario K2E 7L6 Canada
Tel: +1 613-226-8874
Fax: +1 613-226-8639
www.fishersci.ca

Finland

Fisher Scientific Oy
Ratastie 2
01620 Vantaa
Finland
Tel: +358 9 8027 6280
Fax: +358 9 8027 6235
www.fishersci.fi

Ireland

Fisher Scientific Ireland
Suite 3, Plaza 212
Blanchardstown Corporate Park 2
Ballycoolin
Dublin 15
Ireland
Tel: 01 895 5854
Fax: 01 899 1855
www.ie.fishersci.com

Netherlands

Fisher Scientific Emργο B.V.
Zuideinde 70
1121 CM Landsmeer
Netherlands
Tel: +31 20 487 7000
Fax: +31 20 487 7070
www.nl.fishersci.com

Singapore

Fisher Scientific Pte. Ltd.
8 Pandan Crescent
LL4, #05-05 UE Tech Park
128464 Singapore
Tel: +65 6873 6006
Fax: +65 6873 5005
www.fishersci.com.sg

Switzerland

Fisher Scientific AG
Wilstrasse 57, Industrie Wohlen Nord
5610 Wohlen
Switzerland
Tel: +41 56 618 4111
Fax: +41 56 618 4141
www.ch.fishersci.com