MODE D'EMPLOI

CODE ARTICLE	REFERENCE USINE	DESIGNATION
24449	201 425	SONDE à DBO STIRROX G

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Putting into operation / Readiness for measurement

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*Connect probe to the meter.
The probe is immediately ready for measurement.
A polarization of the probe is not required.*

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Note: *StirrOx G operation only possible with connected line adaptor!*

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Recommended ranges of application

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For calibration please refer to operation manual of the instrument.

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For the necessary minimum flow, the probe has been provided with a built-in stirring system.

To turn on the stirring system, press the on/off key.

To turn it off, press the on/off key again.

Recommendation:

After each series of measurements (before storing the probe) clean the probe (exterior cleaning).

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Error symptom	Possible cause	Elimination
Instrument displays 0.0 mg/l or 0 % Sat O ₂ (probe in air)	● No connection between meter and probe ● Defective cable	Control connection between meter and probe Send probe to WTW
Probe cannot be calibrated	● Contaminated exchange cartridge	Refill probe according to operation manual, wait 60 min. and recalibrate
Probe cannot be calibrated even after exchange of electrolyte and membrane head	● Electrode contaminated or probe poisoned	Cleaning of electrode
Leakage is displayed	● Membrane head has not been really tightened ● Hole in the membrane	Screw membrane head tight Exchange membrane head according to operation manual
Wrong temperature indication	● Defective temperature probe	Send probe to WTW
Mechanical damage of the probe		Send probe to WTW
Instrument displays OFL	● Display range exceeded ● Short circuit between gold cathode and lead anode	Regenerate probe; if meter still displays OFL send probe to WTW
Stirrer does not rotate	● Power supply has not been turned on ● Stirring system defective	Turn on power supply Send probe to WTW

Measuring principle	Membrane covered galvanic probe
Temperature compensation	IMT

Measuring conditions


Measuring range	0 - 50 mg/l O ₂
Temperature range	0°C - 50°C
Degree of protection	immersion part of probe: IP 67 whole probe: IP 65
Immersion depth	min. 35 mm max. 70 mm
Operating position	vertical with built-in stirrer
Incident flow	with integrated stirrer

Characteristics when delivered

Zero signal	< 0.1 % of saturation value	
Response time at 20°C	t ₉₀ (90 % of indication of end value after) < 10 s t ₉₅ (95 % of indication of end value after) < 15 s t ₉₉ (99 % of indication of end value after) < 45 s	
Reproducibility	< 0.3 %	
Self consumption at 20°C	0.008 µg h ⁻¹ (mg/l) ⁻¹	
Drift	approx. 3 %/month with permanent polarisation	
Operation time	at least 6 months per electrolyte filling	
Temperature measurement	Temperature adaptation	IMT compensation
	Thermistor type	integrated NTC (30 kΩ / 25°C)
Material	Response time of the thermistor	t ₉₉ < 60 s
	Temperature range	-5 .. +50°C
	Probe accuracy	± 0.2 K
	Membrane head	POM
	Membrane	FEP
	Probe head	PEEK
	Thermistor housing	VA-steel 1.4571
Dimensions	Shaft / housing	POM
	Closing head	Pa 6
	Length of immersion part incl. stirrer	83 mm
	Length of housing	85 mm
	Shaft diameter	43 mm
	Diameter of membrane head	12 mm
	Membrane thickness	13 µm
Weight	Cable length	2 m
	Probe without cable	170 g
	Probe with cable	310 g

Keeping

Storage	in calibration beaker
Temperature range	-5°C ... 50°C

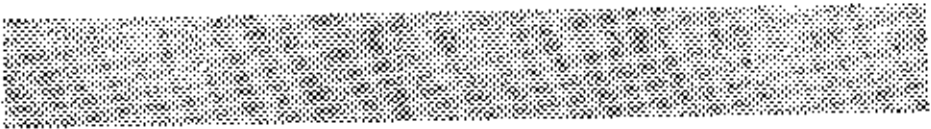StirrOx G