

Thermo Scientific Orion Navigator Pro Measurement Software

System Requirements

Operating System	Microsoft® Windows® 10 (desktop OS) Microsoft Windows 8.1 Microsoft Windows 7
CPU Type	32-bit (x86) or 64-bit (x64) processor @ 1 gigahertz (GHz) or faster
Memory	For 32-bit: 2 GB (Min. 300 MB availability required for the targeted Application usage) For 64-bit: 2 GB (Min. 300 MB availability required for the targeted Application usage)
Visual Reports	Adobe Reader® or equivalent PDF reader program
Display Resolution	1366x768 with True Color
Disk Space	Software Installation 100 MB minimum
Communication Port	One available USB 2.0, USB 3.0 or RS232
.NET Framework	.NET Framework Version 4.5 or higher
Additional Requirements	
Disk Space for application data	1 GB free hard disk Subject to change based on application data growth

Meter Requirements

Compatible Meters

Thermo Scientific Orion Versa Star Pro meters with serial number V10000 or higher and software revision 10.0 or higher (software is user-upgradeable).

Installing Orion Navigator Pro Computer Software

Download the Orion Navigator Pro computer software from the included CD or USB drive.

Note: Refer system requirements table on *page 1*.

1. Double click the **Navigator Pro.msi** file. This will open the Welcome page of the installer.

Name	Date modified	Type	Size
 Orion_Navigator_Pro	02-05-2017 02:54	Windows Installer ...	13,484 KB

Figure 1

2. Click the **Next** button. *Refer to figure 2.*
3. Select the location of the installation folder and then Click the **Next** button.
4. Click the **Next** button to start the installation.

Note: Wait while the software is being installed.

5. Click the **Close** button when the message box confirms that the software has been successfully installed.

6. The Versa Star Pro Navigator icon  should now be available on the desktop. Double click the Versa Star Pro Navigator icon  to launch the application software. *Refer to figure 3.*

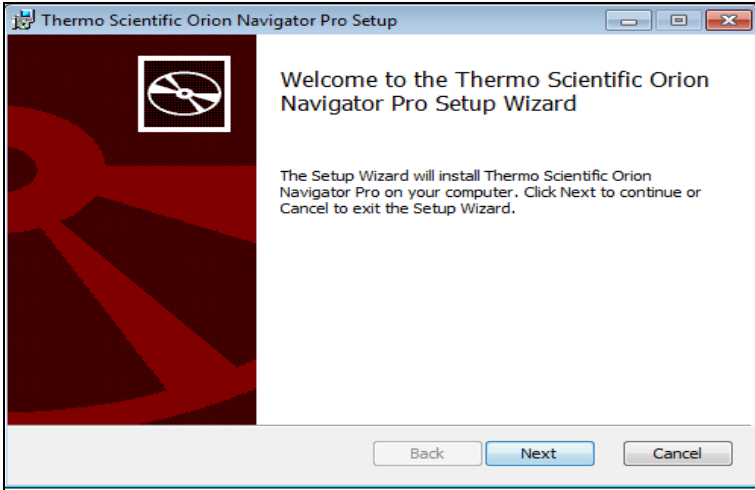


Figure 2

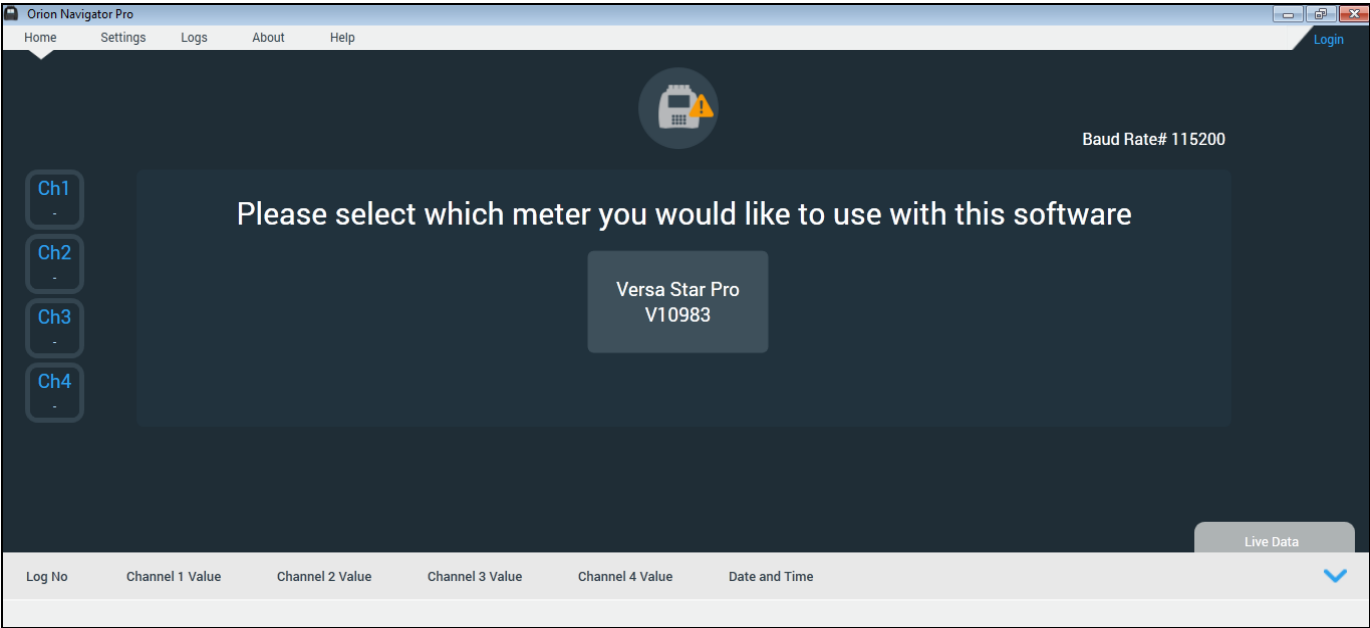


Figure 3

Licensing and Installation

Getting a permanent license key

1. Double click on the application  icon to open the application.
2. Click **Yes** to register the product for the first time.
3. The license registration screen will open. *Refer to figure 4.*
4. Click on the **Copy to Clipboard** tab to copy the UID. *Refer to figure 5.*
5. Click **OK**.
6. Send an email to **SGSI3TS@thermo.com** after pasting the UID to get a license key. *Refer to figure 6.*

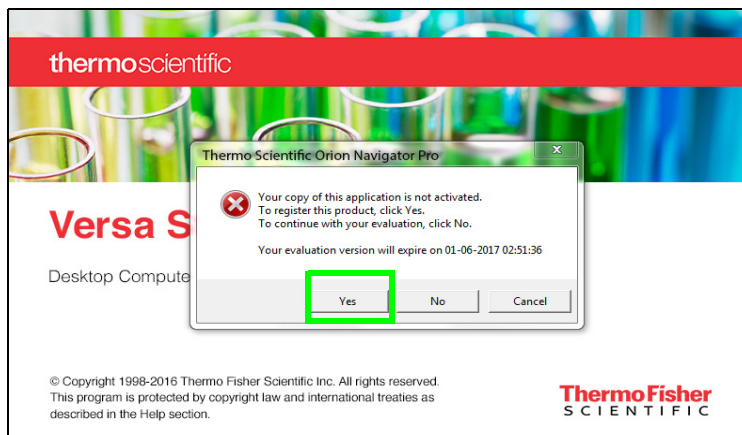


Figure 4

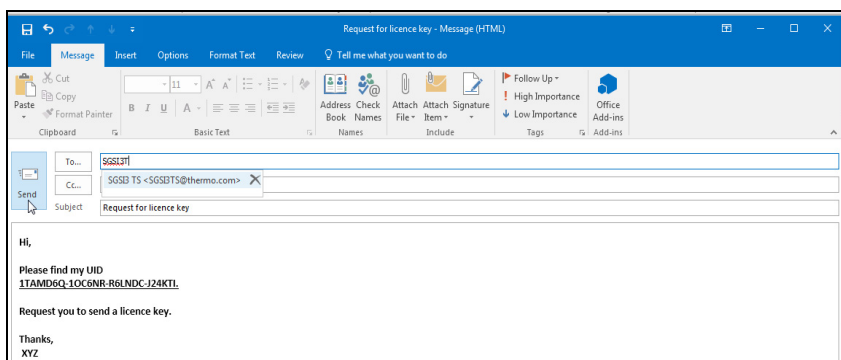


Figure 6

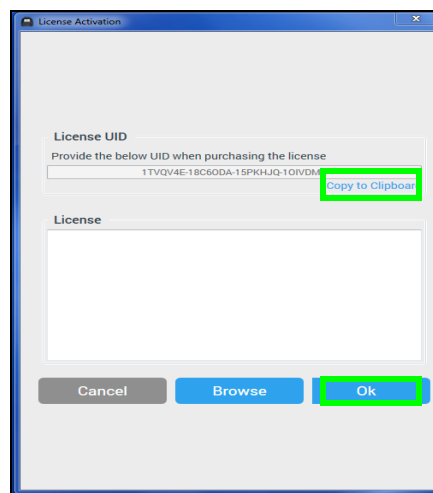


Figure 5

Installing the license

After receiving the license key.

Note: Make sure you save the license key on your system.

7. Create a new folder in your system and place the (.lic) file. *Refer to figure 7.*

Name	Date modified	Type	Size
1YS1HPY-15TV2P8-RZFY47-QK6LDT.lic	31-03-2017 19:59	LIC File	2 KB

Figure 7

8. Right click on the application  icon to **Run as administrator** to install the license key. *Refer to figure 8.*

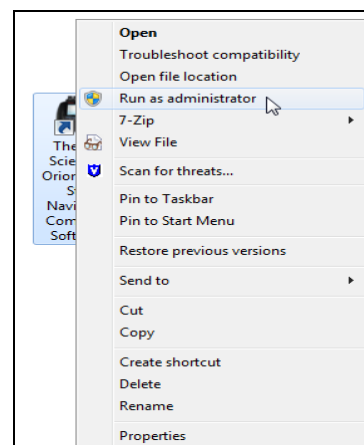


Figure 8

Note: The application will navigate you to the registration screen, *Refer to figure 4.*

9. Click on **Browse** tab to navigate to the license key (.lic) file location. *Refer to figures 9 and 10.*

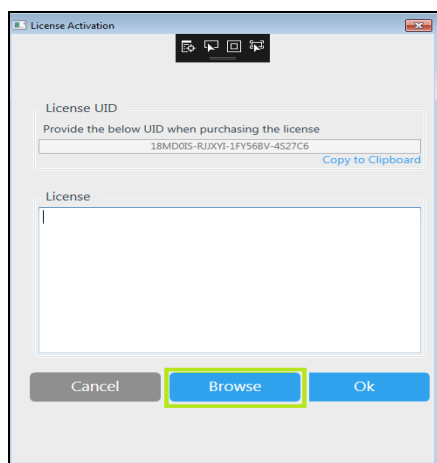


Figure 9

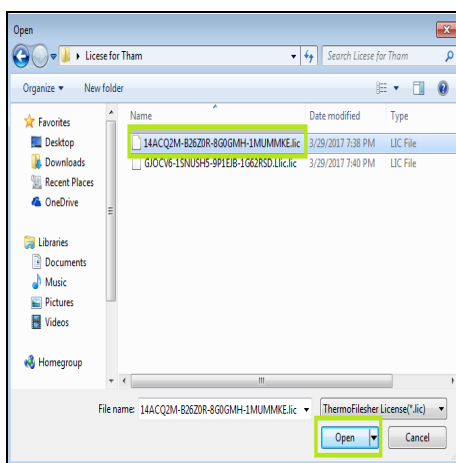


Figure 10

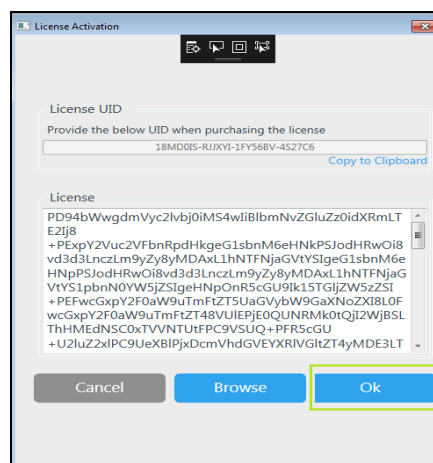


Figure 11

10. Click **OK**. *Refer to figure 11.*

Note: After installation, the application will close automatically.

If you already have a license key

1. Click on **Browse** tab to navigate to the license key (.lic) file location. *Refer to figures 9 and 10.*
2. Click **OK**. *Refer to figure 11.*

thermofisher.com/water

© 2017 Thermo Fisher Scientific Inc. All rights reserved. Microsoft Windows is registered trademark of Microsoft Corporation. Adobe Reader is registered trademark of Adobe Systems Incorporated in the United States and/or other countries. All other trademarks are the property of Thermo Fisher Scientific and its subsidiaries unless otherwise specified.

Water and Lab Products

Australia: (613) 9757-4300 In Australia: (1300) 735-295 **China:** (86) 21-6865-4588
Germany: (49) 6184-90-6321 **India:** (91) 22-6716-2261/2247 **Japan:** (81) 045-453-9175
North America: 1-978-232-6000 Toll Free: 1-800-225-1480 **Singapore:** (65) 6778-6876

Find out more at thermofisher.com/water

68X001440 Rev.A

ThermoFisher
SCIENTIFIC